


# Think You Can't Afford A New Nurse Call System... Think Again.

## *VisionLink™* II Wireless Nurse Call System

### Life Safety Made Easy

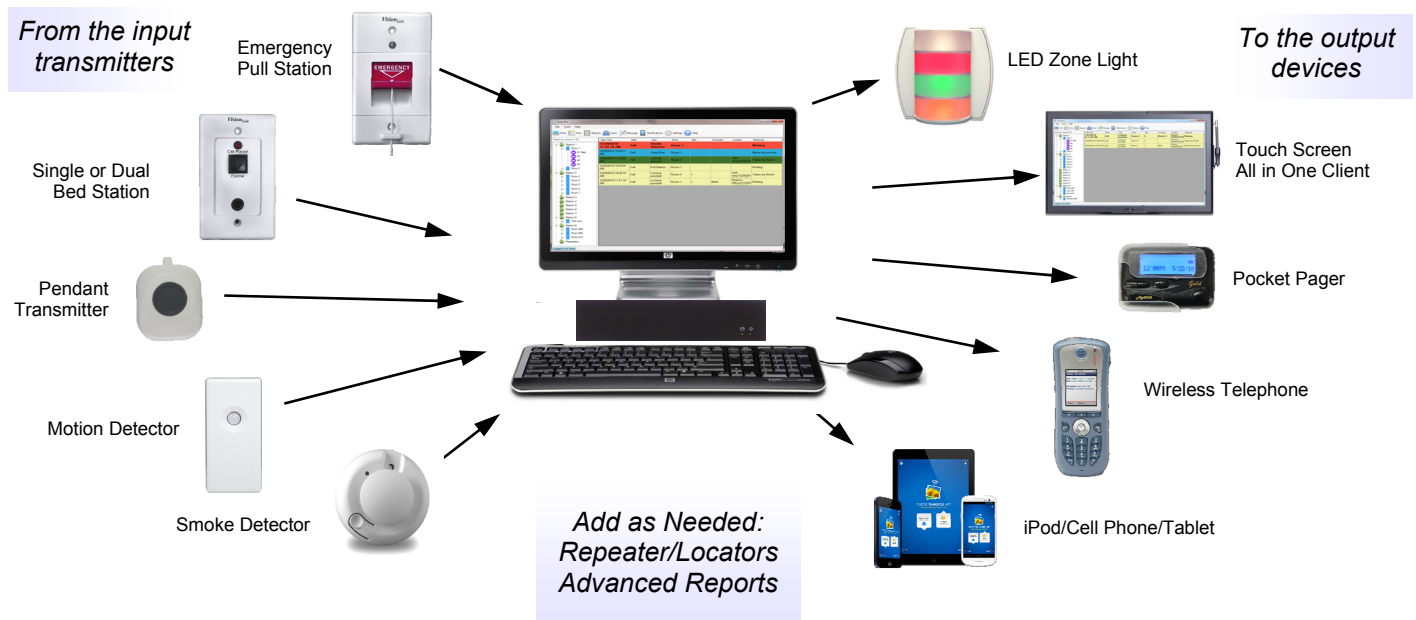

### Start Small, Think Big!

The Vision Link II® Wireless Nurse Call System uses the Windows 10® operating system, 900 MHz spread spectrum wireless technology and an industrial computer with solid state hard drive to link your residents with your staff. Pendant transmitters signal assistance is needed plus provide fall alarm, wander alerts and resident location if desired. The scheduler reminds caregivers to give medications, lock doors or perform other tasks. No expensive licenses or upgrades are required to add more equipment at any time. Wireless repeaters and stations are easily mounted when adding a new wing or building. LED corridor and zone lights are connected on a single pair of wires for fast, easy installation. The Vision Link II® also offers you eight call analysis reports, remote access, email notification and automatic backup at no additional cost.

### Flexible Software

Patient information, statistical call analysis reports with staff response times, active and passive resident and staff check-in, patient bed turn notification, flexible alarm notification with escalation and so much more are all standard features of the powerful and flexible Vision Link II® integration platform. Engineered for the future, new developments are always backwards compatible to ensure it accommodates the growth and requirements of your community.

### Staff Mobility

Caregivers can receive alarm call notification not just on consoles but also on touch screen clients, Android and IOS devices, pocket pagers, two way radios, telephones, email, cell phones, iPod's, iPad's, LED reader boards, remote monitors and more with flexible escalation for unanswered assistance calls. All events are logged including staff response time for password protected management reports. In addition to alarm notification, text messages can be silently sent from any console or computer to any pocket pager, cell phone or employee group eliminating institutional overhead paging. Send staff text messages for calls on hold plus scheduled reminders as needed. Notify staff for resident elopement, door/window exit, smoke detectors, motion detectors, bed/chair exit pads, equipment status, mass emergency notification, whatever your facility requires.

### The 900MHz Advantage

Because the FCC allows higher transmit power for 900MHz spread spectrum systems, Vision Link II provides superior in-building range and up to 10 times the range of 300MHz products. It is also fully supervised and notifies staff when a low battery or inactive device is detected. In contrast, it takes some systems up to 24 hours to report an inoperative station and some wired and wireless nurse call systems are not supervised at all.


# Ensure Safety for Residents and Staff with a *VisionLink™* II Wireless Emergency Call System

## EQUIPMENT OPTIONS

- Location, elopement and fall alarm transmitters
- Wireless patient and code blue stations
- Wireless staff and resident emergency stations
- LED Corridor and Zone Lamps (white or four color)
- Custom sconce indicator lamps
- Alpha-numeric pocket pager & LED reader boards
- Central station alarm monitoring interface
- Mobile App for Android, iPod, iPad, iPhone
- Cisco®, Spectralink®, Ascom & Panasonic interface
- Email, Two way radio, and cell phone interface
- Networked consoles can select response areas
- Remote access for reports and programming
- Door/Window/Bed exit notification
- Wandering patient management
- Smoke/motion/glass breakage detection
- Active or passive resident check-in status
- Repeater/Locators for campus coverage
- Computer, server and desktop consoles
- Locator and Elopement beacons
- Android and IOS App

## Call Today

Our experienced, knowledgeable technical staff would be happy to answer your questions or give you ideas on how to put wireless to work for your facility.

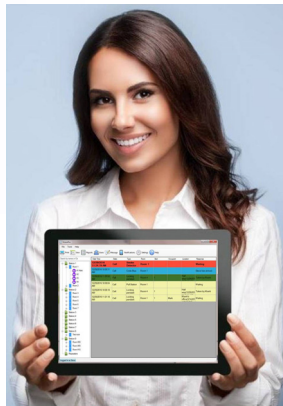

Date Time	State	Type	Room	Bed	Occupant	Location	Response
12/29/2016 11:24:19 AM	Call	Smoke Detector	Room 1				Waiting
12/29/2016 10:58:31 AM	Call	Code Blue	Room 1				Steve has arrived
12/29/2016 11:20:50 AM	Fail	Locking pendant	Room 4	1		Hall way(12/29/2016)	Taken by Khanh
12/29/2016 10:58:04 AM	Call	Pull Station	Room 1				Waiting
12/29/2016 10:58:10 AM	Call	Locking pendant	Room 4	1		Hall way(12/29/2016)	Taken by Khanh
12/29/2016 11:21:10 AM	Call	Locking pendant	Room 1	1	Mark	Khanh's office(01/4/2016)	Waiting

## STANDARD FEATURES

- Industrial computer & solid state hard drive
- 900 MHz spread spectrum technology
- Fully supervised battery and operation
- Alarm calls reset at the source
- Password protection
- Visual/audible alarm prioritization
- User replaceable long-life battery
- Statistical call reports
- Battery back-up included
- Email notification

## MORE FEATURES

- Available ETL Listed to UL1069 Ed. 7 or UL 2560
- Compatible with other call systems
- Pendant fall alarm & wander notification
- Internal & external batteries monitored
- Programmable call escalation
- Advanced management reports
- Silent text messaging
- Quality assurance verification
- Timed egress alarm notification
- Staff presence indication
- Electronic health records interface
- Bed turn scheduler
- Text messaging

Date Time	Room	Bed	Device type	Device ID	Event Type	Responded Person	Response Status	To Room Elapsed Time	In Room Elapsed Time	Take To Attire Elapsed Time	Actions	Note
12/29/2016 7:51:57 AM	Room 4	1	Locking pendant	511	Call	Khanh	Completed	00:29	00:03	00:08	Clean room	
12/29/2016 7:53:56 AM	Room 4	1	Locking pendant	511	Call	Khanh	Completed	01:21	04:25	00:15	Vaccum	
12/29/2016 7:54:55 AM	Room 1		Code Blue	30210	Call	Khanh	Completed	00:27	00:03	00:00	Medicine	
12/29/2016 11:25:53 AM	Room 4	1	Locking pendant	511	Fall	Khanh	Completed	01:29	00:25	01:07	Bath, Medicine	Test 2
12/29/2016 1:03:13 PM	Room 4	1	Locking pendant	511	Call	Khanh	Completed	02:06	00:20	01:14	Meal	
12/29/2016 1:03:16 PM	Room 1		Code Blue	30210	Call	Khanh	Completed	01:52	00:49	00:00	Fall	
12/29/2016 1:03:29 PM	Room 1		Pull Station	30208	Call	Khanh	Completed	01:58	00:44	00:00	Clean room, Fall, Medicine	Test
12/29/2016 1:03:40 PM	Room 1		Smoke Detector	30209	Call	Khanh	Completed	01:54	00:16	00:33	Bed mobility	
12/12/2016 7:32:57 AM	Room 1	2	Locking pendant	325	Call	Mark	Taken	00:00	00:00	00:00		
12/13/2016 12:58:57 PM	Room 4	1	Locking pendant	511	Call	Khanh	Completed	00:17	00:18	00:00	Clean room	



**Systemstechnologies: Transforming Technology Into Solutions**

For more information, call or write:

Toll-Free: 888-826-3394  
 Phone: 208-762-6800  
 Fax: 208-762-4877

Address: 11310 N. Government Way  
 Hayden, ID 83835  
 Email: [sales@wirelessnursecall.com](mailto:sales@wirelessnursecall.com)



Systems Technologies reserves the right to change specifications without notice. Some features may require additional hardware or software. © 2005 Systems Technologies, Inc. All rights reserved. Patents pending worldwide. All products or services mentioned herein are the trademarks of Systems Technologies, Inc., or their respective owners.